

# DynaTrack™ Plus

## patella stabilizer



With its patented design, the DynaTrack Plus aids in the treatment of patellafemoral dysfunction. In addition to pain relief, your patients and athletes will love the comfort and fit the DynaTrack provides.


### Dynamic Buttress System

Unique buttress design applies continuous pressure to the patella during full range of motion of the knee. The straps originating from the buttress allow you to easily adjust the level of tension. The buttress makes direct contact with the patella and is lined with Skinloc™ material to prevent migration.



### Medial and Lateral Plastic Hinges

Removable hinges enhance the brace's ability to conform to the knee during flexion and provide a stable base for the dynamic buttress.

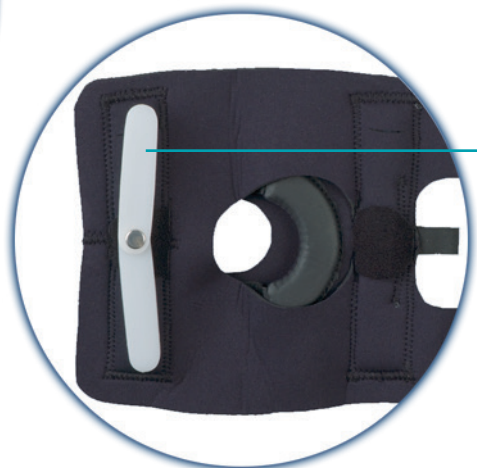

### Open Popliteal

A large popliteal opening permits a high degree of flexion with no discomfort.



### Anterior Closure

The wraparound design allows for ease of application and an optimum fit.

Bilateral design fits left or right knee

Neoprene	Size	Knee Circ.*	CoolFlex
117001	X-Small	10" - 12"	117021
117002	Small	12" - 14"	117022
117004	Medium	14" - 16"	117024
117005	Large	16" - 18"	117025
117006	X-Large	18" - 20"	117026
117007	XX-Large	20" - 22"	117027
117008	3X-Large	over 22"	117028

\*measure circumference of knee at mid patella

CoolFlex version is breathable and lightweight for superior comfort.

Suggested HCPCS: L1820

Medical Specialties, Inc.  
4600-K Lebanon Rd.  
Charlotte, NC 28227  
p: 800-334-4143 f: 800-694-9060  
p: 704-694-2434 f: 704-694-9060  
e-mail: request@medspec.com  
www.medspec.com  
U.S. Patent # 7,004,919

©2006-2017 Medical Specialties, Inc.  
DynaTrack is a trademark of Medical Specialties, Inc.

Neoprene version provides greater compression and warmth.

