


# EVO<sup>®</sup> Ankle Stabilizer



## The next EVolution in ankle bracing technology™

The EVO ankle stabilizer is an EVolutionary step forward in ankle protection. With its patented design, the EVO provides even greater inversion / eversion control and protection against syndesmotic ankle sprains than the original ASO. In addition, the dynamic cuff has been upgraded to offer superior comfort and support.

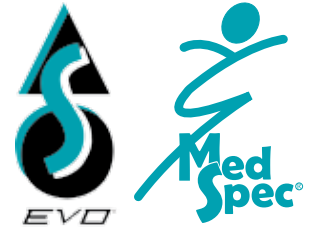

actual product

### Dynamic Cuff

Internal cuff redesigned for lower profile and enhanced circumferential compression around the distal tibiofibular joint. In addition, the dynamic cuff is now lined with CoolFlex™ padding for improved comfort and Skinloc™ material to resist distal migration.

### Stirrup & Stabilizing Straps

Inner stirrup strap is redirected through slots in the dynamic cuff and continues outside the boot to form the stabilizing straps which better capture the calcaneus (heel) and resist inversion or eversion of the ankle to help stabilize the talocrural and subtalar joints.

### CoolFlex Achilles Pad & Tongue

CoolFlex™ padding contains spandex fiber to provide superior comfort along the Achilles tendon and dorsum of foot. The shoe lace is attached to the tongue to ensure lace is always centered.

### Invisible Seam

For improved comfort and durability, there is no seam under the foot.



A portion of the elastic cuff has been removed and white straps used in the above photo for illustrative purposes.

Bilateral design fits left or right foot

Black	White	Size	Ankle Circum.*
264091	264081	x-small	10" - 11"
264092	264082	small	11" - 12"
264094	264084	medium	12" - 13"
264095	264085	large	13" - 14"
264096	264086	x-large	14" - 15"
264097	264087	xx-large	15" - 16"

\*measure ankle circumference according to diagram below



PDAC Verified HCPCS:  
L1902

### Indications:

- Inversion or eversion ankle sprains
- Syndesmotic ankle sprains (high ankle sprains)
- Medial or lateral ankle instability
- Subtalar joint instability

Men's Shoe Size	6-7	7 1/2-9	9 1/2-11	11 1/2-13	13 1/2-15
Narrow Foot		xs	s	m	l
Regular Foot	xs	s	m	l	xl
Wide Foot	s	m	l	xl	xxl
Women's Shoe Size	5-6	6 1/2-8	8 1/2-10	10 1/2-12	12 1/2-14
Narrow Foot		xs	s	m	l
Regular Foot	xs	s	m	m	l
Wide Foot	s	m	l	l	xl

Medical Specialties, Inc.

4600-K Lebanon Rd. Charlotte, NC 28227 p: 800-334-4143 f: 800-694-9060 p: 704-694-2434 f: 704-694-9060

e-mail: request@medspec.com www.medspec.com

Patent Numbers: US 8,721,578; US 8,808,215; US 9,375,339; US 9,737,430; US 10,383,404 AU 2008317265;

CA 2,703,266; European (AT, BE, CH, DE, DK, FR, GB, IE, NL, NO, SE) 2203135; patent pending

ASO and EVO are registered trademarks of Medical Specialties, Inc.

© 2008-2023 Medical Specialties, Inc.



rev 5/23